

Installation Guide


Mounting Kit for Mounting Philips Flexible Module Server (FMS) or Docking Station on 2" Post

The purpose of this guide is to describe mounting of FMS or docking station on 2" post.

Parts Reference

The following parts and hardware are included with this installation kit (hardware not shown):


Item #	Description	Qty
1	Mounting Adapter with Front Clamping Block	1
2	Clamping Block, Rear Half	1
3	M4 x 12 mm Pan Head Machine Screw (PHMS)	2
4	M4 x 16 mm PHMS	2
5	#10-32 x 9/16" PHMS	2

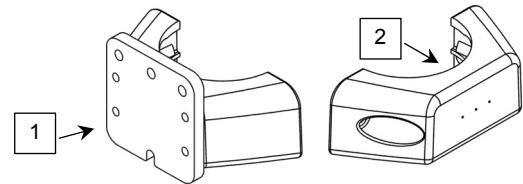


Tools Required

- Phillips screwdriver (not provided)


Attaching Mounting Adapter to FMS

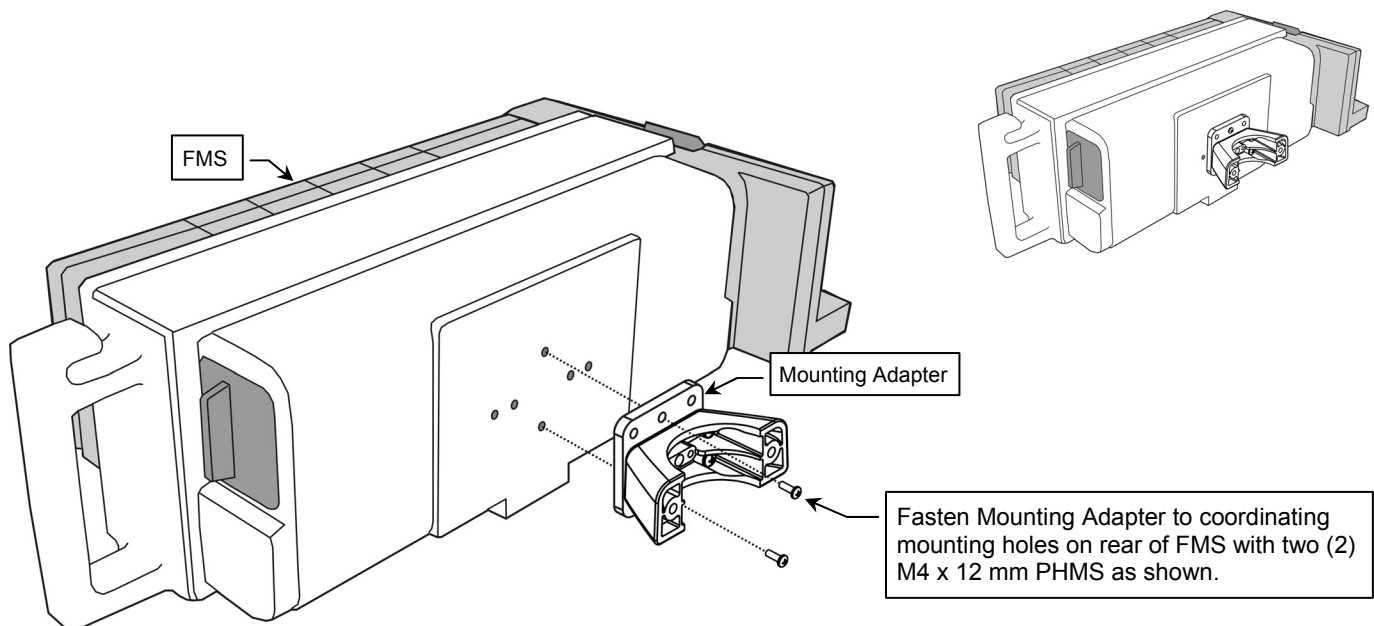

Attaching Mounting Adapter to Docking Station

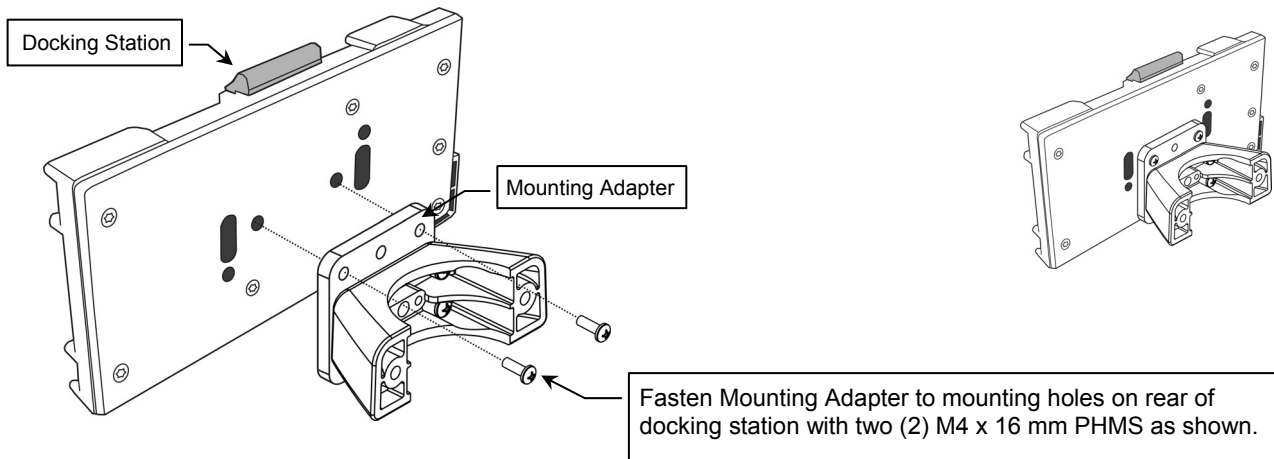

Mounting FMS or Docking Station on 2" Post

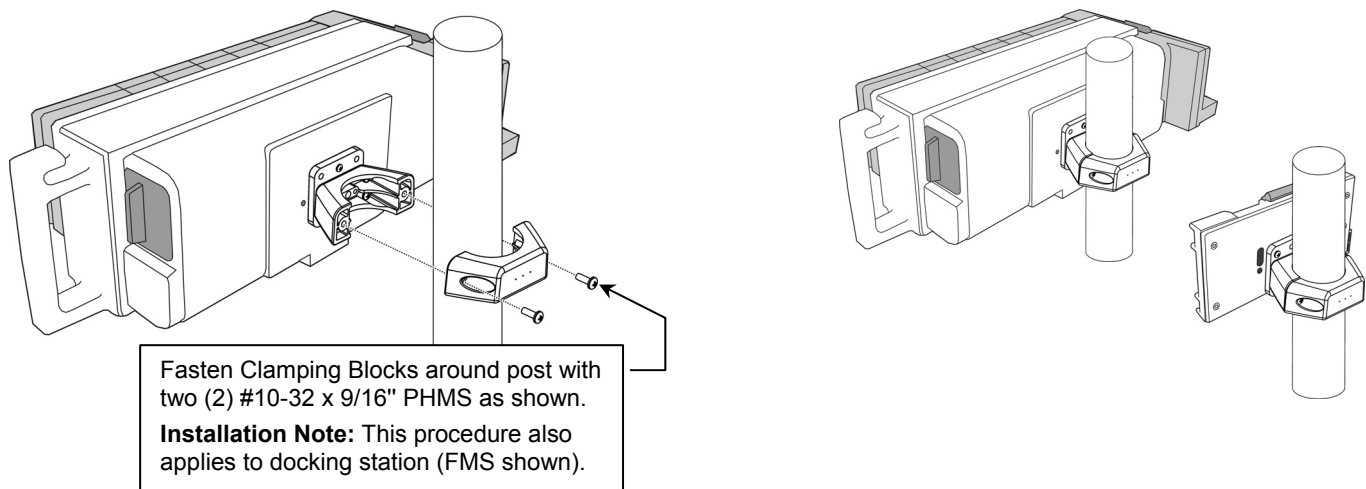

Routine Maintenance

Periodically inspect all mounting hardware. Tighten or adjust as necessary for optimal operation and safety.

Cleaning the Mounting Assembly

CAUTION: GCX makes no claims regarding the efficacy of the listed chemicals or processes as a means for controlling infection. Consult your hospital's infection control officer or epidemiologist. To clean or sterilize mounted instruments or accessory equipment, refer to the specific instructions delivered with those products.

1. The mounting assembly may be cleaned with most mild, non-abrasive solutions commonly used in the hospital environment (e.g. diluted bleach, ammonia, or alcohol solutions).
2. The surface finish will be permanently damaged by strong chemicals and solvents such as acetone and trichloroethylene.
3. Do not use steel wool or other abrasive material to clean the mounting assembly.
4. Damage caused by the use of unapproved substances or processes will not be covered by warranty. We recommend testing of any cleaning solution on a small area of the mounting assembly that is not visible to verify compatibility.
5. Never submerge or allow liquids to enter the mounting assembly. Wipe any cleaning agents off the mounting assembly immediately using a water-dampened cloth. Dry all mounting assemblies thoroughly after cleaning.